

## Anti-Flag Tag Mouse Monoclonal Antibody

**MESGEN**  
INNOVATION BIOTECHNOLOGY

**Catalog Number :** MAN2001

**Lot Number :** Refer to vial

**Packaging Size :** 100uL

### Background

Recognizes Fusion proteins from all species containing the DYKDDDDK epitope

### Overview

**Description :** Mouse monoclonal antibody to Flag-tag [3A6]-Epitope Tag Antibody

**Reactivity :** All

**Tested applications :** WB, ICC/IF, IP

**Formulation :** Liquid


### Properties

**Immunogen :** Synthetic peptide DYKDDDDK conjugated to KLH.

**Clonality :** 3A6

**Isotype :** Mouse IgG1

### Application notes

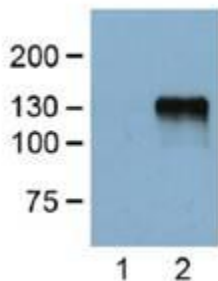

**Fig.** WB analysis (1:5000)

Lane 1: Untransfected control.

Lane 2: HEK293 cells transfected with DYKDDDDK-tagged protein vector at 1:5000.

WB(1:1000-1:10000), ICC/IF (1:100-1:500), IP (1:50-1:200).

Not yet tested in other applications. Optimal dilutions should be determined by the end user.

### Storage instruction

Store at +4°C short term (1-2 weeks). Aliquot and store at -20°C or -80°C

### Other package

Package	List Price (RMB)	Discount Price (RMB)
200ul	2205.00	
500ul	5245.00	
1ml	9802.00	

### Distributor in China ( Mainland )

**Company :** Shanghai hongsheng Biotechnology Co.,Ltd

**Address :** Room 2-205, Xingzhan building, No.185, Liuying Road, Zhabei District, Shanghai

**Tel :** 86-21-56620378 **Fax :** 86-21-56620378

**Email :** [sales@hongshengbio.com](mailto:sales@hongshengbio.com)

**Website :** [www.hongshengbio.com](http://www.hongshengbio.com)

**Note:** The product listed herein is for research use only and is not intended for use in human or clinical diagnosis. Suggested applications of our products are not recommendations to use our products in violation of any patent or as a license. We cannot be responsible for patent infringements or other violations that may occur with the use of this product.

