

Anti-6X His Mouse Monoclonal Antibody

Technical literature is available at: www.mesgenbio.com. E-mail MesGen Technical Services if you have questions on use of this system: tech@mesgenbio.com



Catalog Number : MAN2004

Lot Number : Refer to vial

Packaging Size : 100µL □ 200µL □ 500µL □ 1000µL □

Background

This antibody recognizes all proteins which contain either an N or C-terminal HHHHHHGS epitope tag. It has been shown to recognise HHHHHH alone.

Overview

Description : Mouse monoclonal antibody to 6X His-Tag [2B6]

Tested applications : WB(1:500-1:5000) ICC/IF(1:50-1:500) IP(1:50-1:200)

Formulation : Liquid ,100 µg at 1mg/ml

Properties

Immunogen : 6x His synthetic peptide

Clonality : 2B6

Isotype : Mouse IgG2a

Application notes

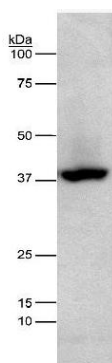


Fig. Lane 1: lysates prepared from HEK-293 cells transfected with his-tagged protein at 1 : 1000 dilution

Optimal working dilutions should be determined experimentally by the investigator.

Storage instruction

Store at +4°C short term (1-2 weeks).

Aliquot and store at -20°C or -80°C

Note: The product listed herein is for research use only and is not intended for use in human or clinical diagnosis. Suggested applications of our products are not recommendations to use our products in violation of any patent or as a license. We cannot be responsible for patent infringements or other violations that may occur with the use of this product.

Distributor in China (Mainland)

Company : Shanghai hongsheng Biotechnology Co.,ltd

Tel : 86-21-56620378 **Fax :** 86-21-56620378

Email : sales@hongshengbio.com

Website : www.hongshengbio.com