

# Goat Anti- Rabbit IgG (H+L)-PE Conjugated

Technical literature is available at: [www.mesgenbio.com](http://www.mesgenbio.com). E-mail MesGen Technical Services if you have questions on use of this system: [tech@mesgenbio.com](mailto:tech@mesgenbio.com)

**MESGEN**  
INNOVATION BIOTECHNOLOGY

**Catalog Number** : MAN4524

**Lot Number** : Refer to vial

## Background

The antibody reacts with whole molecule Rabbit IgG. No antibody was detected against non-immunoglobulin serum proteins.

## Conjugation

PE. Ex: 488nm, Em: 575nm

## Overview

**Description** : Goat polyclonal Secondary Antibody to Rabbit IgG - H&L (PE)

**Reactivity** : Rabbit

**Formulation** : Liquid

**Tested applications** : FC (Flow Cytometry), IF (Immunofluorescent)

IF: 1/50 ~1/1000

Since applications vary, the appropriate dilutions must be used at an assay dependent concentration.

## Properties

**Immunogen** : Rabbit IgG whole molecule

**Clonality** : Goat Polyclonal


**Isotype** : Goat Polyclonal IgG

## Storage instruction

Store at +4°C. and protected from prolonged exposure to light. Do not freeze.

**Note:** The product listed herein is for research use only and is not intended for use in human or clinical diagnosis. Suggested applications of our products are not recommendations to use our products in violation of any patent or as a license. We cannot be responsible for patent infringements or other violations that may occur with the use of this product.

## Application Images

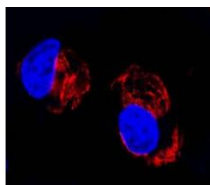

Immunofluorescence Staining of HepG2 cell with primary antibody beta-tubulin and PE labelled goat anti-rabbit IgG as secondary antibody diluted at 1/200

## Distributor in China ( Mainland )

**Company** : Shanghai hongsheng Biotechnology Co.,ltd

**Tel** : 86-21-56620378    **Fax** : 86-21-56620378

**Email** : [sales@hongshengbio.com](mailto:sales@hongshengbio.com)

**Website** : [www.hongshengbio.com](http://www.hongshengbio.com)

MesGen Biotechnology

[www.mesgenbio.com](http://www.mesgenbio.com)